

Wexis Cold Plate CP380



Wexis CP380 maintains an ideal cooling temperature for the tissue blocks pending for sectioning. It can also work to support an embedding machine; or with a wax dispenser to become a simple embedding system.

Design Patent Number: ZL 2014 3 0316960.X
Utility Model Patent Number: ZL 2014 2 0442375.9

Feature

- Eco-Compressor refrigeration system, long service life & efficient freezing;
- Effectively maintaining paraffin block temperature for smooth sectioning;
- Micro-processor temperature control, easy to operate; & superb performance;
- Cold Plate is made of quality aluminum alloy, good conductivity ~ temperature performance;
- Cold Plate with unique elevated edges to avoid water drain off, ensuring clean work surface;
- Energy saving, after reaching the set temperature, system will be paused, saving power and readily stand-by to support daily routine.

Wexis


*Image is for illustration purposes only.

Wexis Group Limited is a global pathology supplies and medical equipment manufacturer. We have advanced professional equipment, dedicated research development laboratory, and production workshop, and are one of the leading manufacturers of pathology supplies in China.

WX-CY-Y051401E

Wexis Cold Plate CP380

SMALL INSTRUMENT SERIES

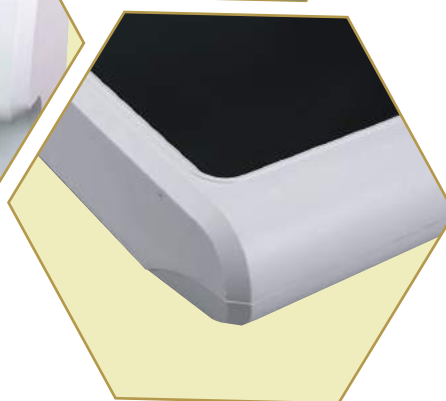

Spacious Cold Plate holds a max. of 48 paraffin blocks at one time;


Top of the instrument provides a platform for placing small accessories & tools.

Unique outer casing with ergonomic designed, no sharp corners, safe and comfortable operation;

SPECIFICATIONS

Physical Size:	535*330*358mm
Cold Plate Area:	270*250mm
Net Weight:	20Kg
Line-in (Voltage):	230V ~ 50HZ or 110V ~ 60Hz
Rating:	170W
Temperature Range:	-2°C~-15°C
Time to Reach the Set Temperature:	≤15min

ORDERING DETAILS

Order No.:	380-00-00
Model No.:	CP380
Description:	Cold Plate CP380
Includes	Main Unit (1) Power Cable (1) Instruction Manual (1)
Certification	Quality Certified Organization CE 

WEXIS SMALL INSTRUMENT SERIES



Floatation/Drying Workstation FL220



Slide Drying Hotplate DR230



Wax Trimmer TR302



Wax Dispenser WD360

Wexis Group Limited

Address: Unit1611, 16/FI, Yuen Long Trading Centre,
33 Wang Yip Street, Yuen Long, NT, Hong Kong
Tel: 852-2567 2882 Fax: 852-2567 6333

Guangzhou Wexis Biotech Limited

Address: No.4 of Yin Sha Road 2, Sha Pu Yin Sha Industrial Park,
Xin Tang, Zeng Cheng District, Guangzhou, P.R, China
Zip-Code: 511388
Tel: 8620-8779 3718 Fax: 8620-8779 3728

Website: www.wexisgp.com
Enquiry: sales@wexisgp.com

*Image is for illustration purposes only.

