
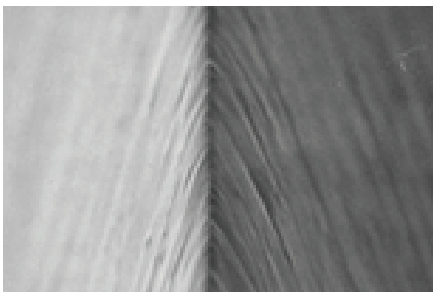


Wexis Microtome Blade

Wexis

The latest Wexis microtome blade is made of the finest carbon steel material. Produced through professional manufacturing process with the latest technologies. These latest models have received numerous appraisals in its sharpness and durability.

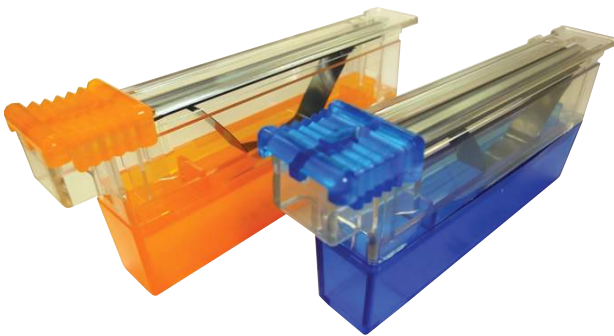

Observation on blade edge [SEM x20,000]



Easy to cut large amount of paraffin blocks

The blade edges are treated by ion cleaning and coated with Teflon (PTFE) that the blades have incredible sharpness and durability. One blade can easily cut 50 to 60 biopsy blocks or 35 to 40 routine blocks in 3 μ m cutting sections.

Cutting hard tissue becomes much easier now SS35 can easily cut hard tissues like uterus, prostate, mammary, thyroid glands, etc

Easy to cut with original form

DS35 can cut sections as the same size as the blocks, eliminate transformations.

Easy to put on slide glass

Special treatment prevents curling of specimen so that a section is easily picked up by the slide glass.

Order Number

Order Number	6011351	6012351
Name	Microtome Blade - DS35 Dura Sharp	Microtome Blade - SS35 Super Sharp
Packaging	50pcs/case,10case/box (500pcs/box)	50pcs/case,10case/box (500pcs/box)
Specification	35°/80X8X0.254mm	35°/80X8X0.254mm



DS35 is suitable for all routine tissue pathology sectioning, fits on all different microtome brands and models, offers perfect durability and good cutting performance.



SS35 offers ideal sharpness, especially suitable for hard tissue sectioning; with good durability.

*Image is for illustration purposes only.

Company Information

Wexis Group Limited is a global pathology supplies and medical equipment manufacturer. We have advanced professional equipment, dedicated research development laboratory, and production workshop, and are one of the leading manufacturers of pathology supplies in China.

We've accumulated practical experience in sales and production for decades. Combining the latest and advanced technologies in the global market with modern management system, we offer to our customers the best quality service, and simultaneously maintaining highly efficient, stable and reliable production.

Wexis Group Limited
Address: Unit1611, 16/FI, Yuen Long Trading Centre,
33 Wang Yip Street, Yuen Long, NT, Hong Kong
Tel: 852-2567 2882 Tel: 852-2567 6333

Guangzhou Wexis Biotech Limited
Address :No.4 of Yin Sha Road 2, Sha Pu Yin Sha Industrial Park,
Xin Tang, Zeng Cheng District, Guangzhou, P. R, China
Zip-Code: 511388 Tel: 8620-8779 3718 Fax: 8620-8779 3728

www.wexisgp.com

Wexis